

# 18

## Shower trolleys



## CHINESPORT TOPICS

Professional tables	1
Electro-medical equipment	2
Assistive and passive training	3
Active exercises	4
Proprioception	5
Pulley therapy	6
Occupational therapy	7
Standing and mobility	8
Tilt tables	9
Parallel bars and staircases	10
Walking aids	11
Treadmills	12
Tractions	13
Medical and postural gymnastics	14
Hydrotherapy	15
Posture analysis	16
Hoisting systems	17
<b>Hygiene solutions</b>	<b>18</b>
Lifting armchairs	19
Patient transport	20

Edition No. 3 - July 2020 - © Copyright Chinesport spa - Italy - Since 1976

At our sole discretion, the company reserves the right to change dimensions and type of construction, and to make improvements and other changes to its products. All reproduction rights of all or part of the designs and illustrations are reserved worldwide. The printing process can't give a perfect reproduction of the colours.

Chinesport thanks all those who contribute to the development of the contents of this document.



UNI EN ISO 9001:2015  
UNI EN ISO 13485:2016

# Shower trolleys



<b>1. Introduction</b>	<b>04</b>
<b>2. Configuration options</b>	<b>05</b>
2.1 Type of control	06
2.2 Inside width	07
2.3 Length extension	08
2.4 The backrest	08
2.5 Wheels and braking system	09
<b>3. Accessories</b>	<b>10</b>
3.1 New shower mattress	10
3.2 Posture cushions	11
3.3 Electric parts	12
3.4 Other accessories	12
<b>4. Models</b>	<b>14</b>
4.1 LD31 for paediatric use	14
4.2 LD32 for adults	15
4.3 LD33 for bariatric use	16
<b>5. Showering in bed</b>	<b>17</b>
<b>6. Appendixes</b>	<b>18</b>
6.1 Spare parts list	18
6.2 Exploded views	21

# Shower trolleys *Introduction*

All our shower trolleys have standard additional function to also adjust them to Trendelenburg position.



Our shower trolley are developed in collaboration with:




**Centre hospitalier  
universitaire vaudois**  
Lausanne - Switzerland

## THE QUALITY THAT YOU CAN SEE ON THE LONG TERM



SIGNS OF CORROSION SAMPLE

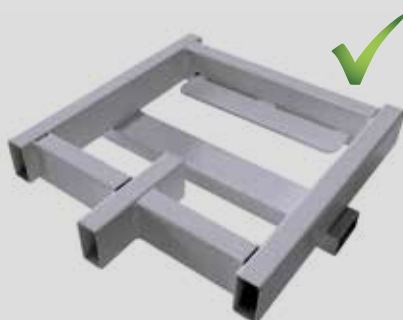

OUR TREATED SAMPLE

The picture beside here shows one of our samples previously treated with special anti-corrosion treatment, after a long test with exposure to a critical environment in a specialised laboratory. The sample is compared to an untreated piece which, after the same tests, shows significant damage due to corrosion agents.



Our shower and transfer trolleys for children, adults and bariatric patients, are suitable for use both in care facilities and at home. We have devised solutions that make it possible to adapt the shower trolley according to the preferences of caregivers or family members and to the user's anthropometric characteristics.


## For paediatric use



**L D 3 1**

SHOWER TROLLEY 31

The product code for Chinesport shower trolley is an "open code" and each code position has a particular meaning, which corresponds to one or more features of the shower trolley.

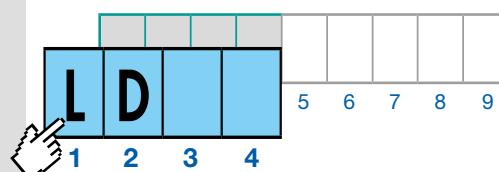

### CODE POSITIONS 3-4

#### The basic models

These two first positions of the code identify the different basic models such as for children, adults or bariatric use.

### CODE POSITION 5

#### Height adjustment

- C** - Hydraulic, with double pedal
- D** - Electric, with hand control
- G** - Electric, footswitch with its holder

### CODE POSITION 6

#### Inside width

- 1** - Inside width 60 cm
- 2** - Inside width 68 cm
- 3** - Inside width 80 cm

### CODE POSITION 7

#### Extension length

- 1** - Fixed rails on head and foot end
- 2** - Head end rail partially folding
- 3** - Partially folding rails on head and foot end

### CODE POSITION 8

#### Tilt adjustable backrest

- 1** - Non-Tilting backrest
- 2** - Tilting backrest

### CODE POSITION 9

#### Wheel system for transfer

- 1** - Wheels with individual brake
- 2** - Wheels with central braking system

## For adults



**L D 3 2**

SHOWER TROLLEY 32

## For bariatric use




**L D 3 3**

SHOWER TROLLEY 33

# Shower trolleys *Configuration Guide*

## POSITIONS 5

### Height adjustment and type of control

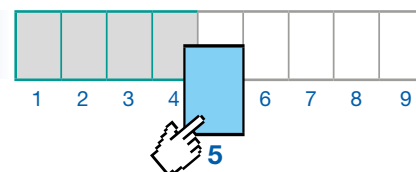

The shower trolley height can be adjusted via electric actuator or hydraulic piston. The electric version can be adjusted using the pedal control or handset. Electric versions are supplied with battery charger.



- C Hydraulic, with double pedal**  
Adjustment via hydraulic piston, with pedal control available on both sides of the shower trolley. This hydraulic adjustment is ideal for situations where there is no electric power supply or is preferred for minimum periodical maintenance.



- D Electric, with hand control**  
Ergonomic hand control with magnet. Caregivers can attach it wherever they wish on the frame metal parts. It comes with connection cable.




- G Electric, with foot switch on its holder**  
The control is securely placed in a special holder, which has to be attached at the base of the structure. The pedal control on its holder makes the transferring of the shower trolley easier.



We present the three basic models for pediatric, adult and bariatric use with an overview of the various options available. The video shows the capabilities of the Chinesport shower trolleys with details and enhance the configurability of the product.



## POSITIONS 6

### Inside width of the shower mattress


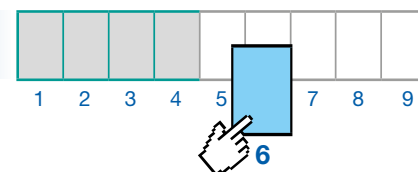


12



Various mattress inside widths are available to fit the end user's anthropometric characteristics and the consequent needed space for a comfortable use. In particular there are three standard inside width options for the mattress: for paediatric use, for adults or bariatric use. In addition it's possible to choose option of partial foldable siderails at the head and foot side.

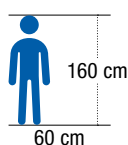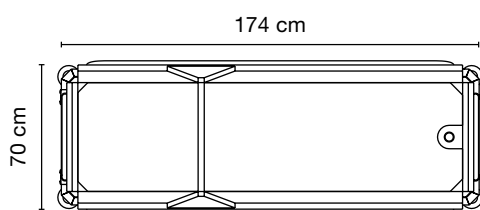

- 1 60 cm inside width mattress**  
The overall dimensions of the shower mattress are 174 cm in length and 70 cm in width. Consequently the inside dimensions are 160 cm length and 60 cm width. That's the available surface for the end user. This mattress option is suitable in particular way for children up to 160 cm tall. Weight: 4 kg

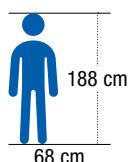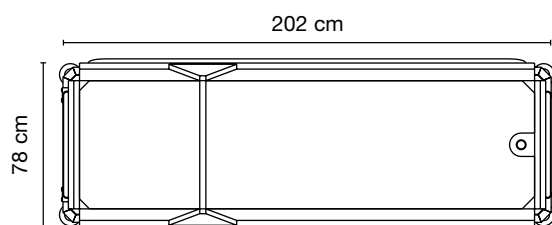

- 2 68 cm inside width mattress**  
The overall dimensions of the shower mattress are 202 cm in length and 78 cm in width. Consequently the inside dimensions are 188 cm length and 68 cm width. That's the available surface for the end user. This mattress option is suitable in particular way for adults up to 188 cm tall. Weight: 4.5 kg

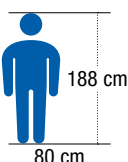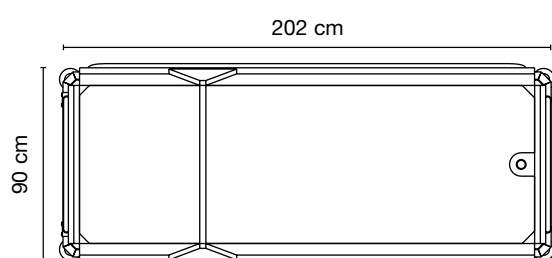


- 3 80 cm inside width mattress**  
The overall dimensions of the shower mattress are 202 cm in length and 90 cm in width. Consequently the inside dimensions are 188 cm length and 80 cm width. That's the available surface for the end user. This mattress option is suitable in particular way for obese or large people and up to 188 cm tall. Weight: 5 kg

# Shower trolleys *Configuration Guide*

## POSITIONS 7

### Mattress length extension

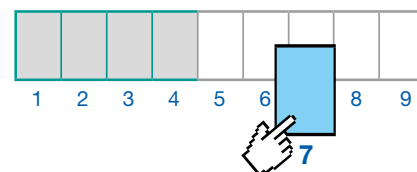


The length of the mattress usable space can be extended by tilting the head end rail partially, and tilting the foot end rail as well if required. The carer can then easily lift and tilt the siderail up to a 40° angle as needed. The result is to gain 7 or 14 cm in length to adapt the shower mattress to people taller than 188 cm.



- 1 Head and foot safety side rails not foldable**  
Safety siderails on head and foot end can not be folded. The mattress inside length remains 188 cm.
- 2 Head safety siderail foldable**  
The safety siderail on head side is partially foldable of 40° to extend the mattress length of 7 cm. So in total the mattress length is 195 cm for a tall end user.
- 3 Both head and foot safety side rails foldable**  
Both safety siderails on the head and foot end can be partially folded of 40°. The mattress inside space is extended in such a way up to 14 cm and in total the inside length becomes of 202 cm.

## POSITIONS 8

### Tilt-adjustable backrest

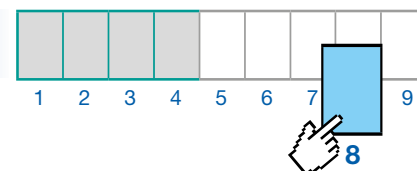


The shower trolley backrest can be tilt adjustable to give greater comfort to both end user and caregiver during washing, as it allows to change the patient's position. The backrest can be easily tilted by means of gas spring with ergonomic lever.

- 1 Non-tilting backrest**  
The patient remains in a supine position inside the mattress. The patient's head can be lifted by using a cushion that is an accessory to be ordered separately.
- 2 Tilting backrest**  
The backrest can be tilt adjustable in positive up to 63°. The tilt can be locked or modified in any intermediate position. A separate extension is an accessory which can also be requested to provide a further head support.



## POSITIONS 9

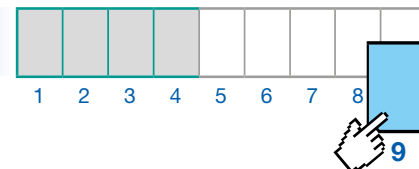

12

### Wheels and braking system

You can choose between individually braked wheels or central braking system via a pedal. In both cases there is also a “directional lock” system for easier transportation on long, straight stretches. All Chinesport shower trolleys are fitted with swivelling bumpers integrated into the deck structure.



Wheel with individual brake



Wheel for directional lock

- 1 Wheels with individual brake**  
The wheels have 12.5 cm diameter and all provided with individual brake. In addition one of them is identified with green part that enables the directional function.



Activated central brake




Activated directional lock

- 2 Wheels with central braking system**  
Wheels have 15 cm diameter with central braking system and “directional lock” function. These functions can be activated alternately using the same pedal system, that's available on both short sides of the trolley.

# Shower trolleys *The shower mattress*

## A new light and soft mattress

As result of a long research and testing phase our Research & Development department has released in production a new standard mattress for all three sizes of our shower trolleys: for pediatric, adults or bariatric use. The new mattress is in a soft material, smooth and elastic at the touch, and allows to make faster and efficient the cleaning and disinfection operations for its light weight.

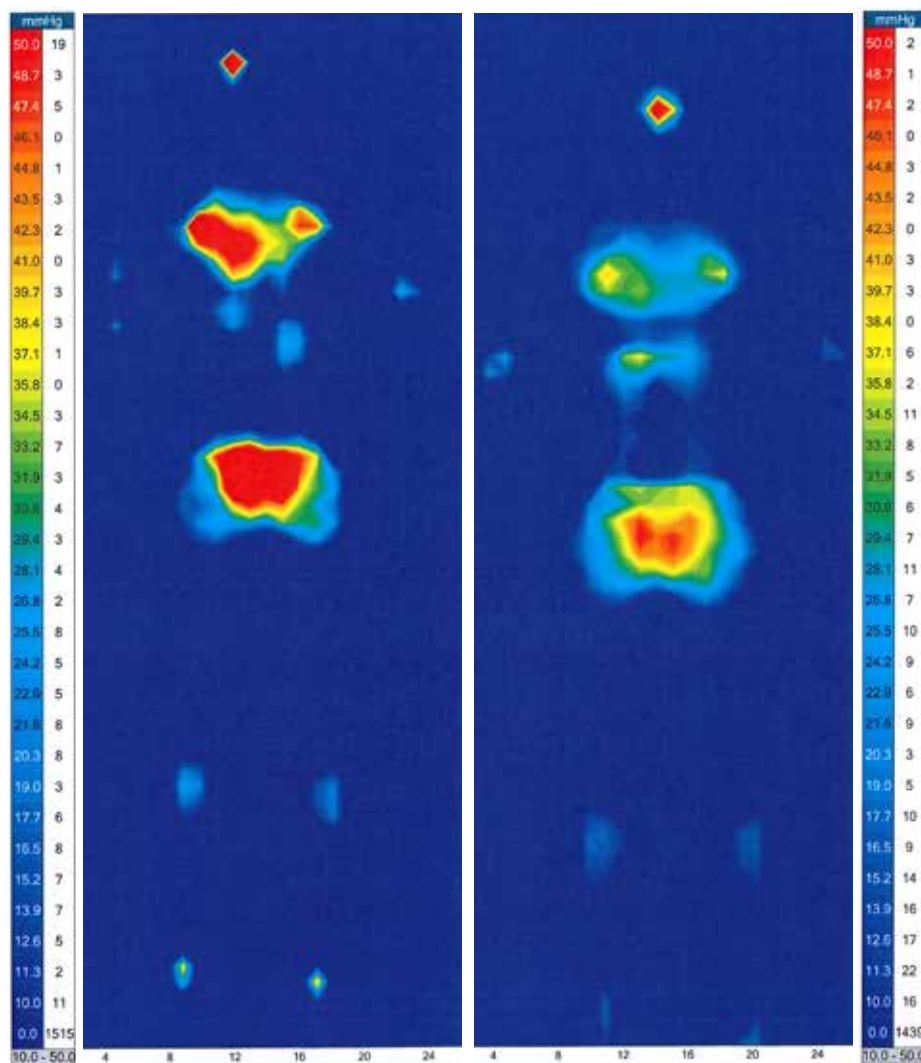

The change of material also enables more comfort for the patient, having introduced a different thickness for both the bottom of mattress and the rising edges. Laboratory test by using sensors have allowed to realize a significant reduction of the medium pressure, the peak of pression and a wider body surface in contact with the mattress compared with other classic ones.

PARAMETERS OBSERVED	TRADITIONAL MATTRESS	NEW GENERATION MATTRESS
Medium pressure [mmHG]	31.5	23.1
Peak of pressure [mmHG]	113.6	69.7
Body surface contact [in²]	232.81	351.56



### Ergonomic controls

Our shower trolley have ergonomic controls to adjust the head end tilt, and to adjust the platform tilt angle for water drainage as well as Trendelenburg position. All these adjustments are assisted by reliable gas springs.

## SHOWER MATTRESS ACCESSORIES

12



To whom are already using our shower trolleys and in case of need to replace the existing mattress, they can buy as spare part the mattress with the new features.

### SP0041 SOFT MATTRESS - SMALL

Mattress inside dimensions are 60 cm in width and 160 cm in length. This mattress can be requested for model LD31. Weight: 4 kg

### SP0035 SOFT MATTRESS - MEDIUM

Mattress inside dimensions are 68 cm in width and 188 cm in length. This mattress can be requested for models LD32 and LD33. Weight: 4,5 kg

### SP0042 SOFT MATTRESS - LARGE

Mattress inside dimensions are 80 cm in width and 188 cm in length. This mattress can be requested for models LD32 and LD33. Weight: 5 kg



### AC1237 SEMI-CYLINDRICAL CUSHION

Dimensions: 55 x 15 x 7 cm



### AC1236 WEDGE CUSHION

Dimensions: 55 x 55 x 14/2 h cm.  
15° inclination.



### AC0574 BACKREST EXTENSION

For greater support to the backrest. It is only available for the versions with inside width 68 cm and reclining backrest option.


### AC1126 SHOWER MATTRESS STRAPS

Set of straps for quick application to all four sides of the shower mattress. These straps guarantee further stability of the raised edges of the mattress to be fixed to the side rails. The set is made of 2 short straps with clips for long sides and 2 long straps for short sides.






## ELECTRICAL ACCESSORIES

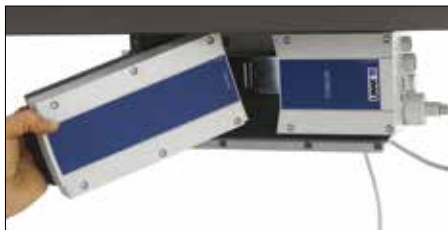

### AC0071 SPARE BATTERY

This is a rechargeable, 24 V spare battery recommended to prevent any interruptions during use.

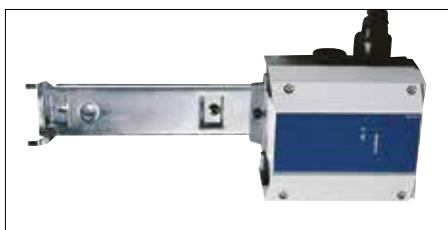

### AC0318 WALL BATTERY CHARGER

This is to recharge both the main and spare battery by connecting it to the mains power.



### AC0578 SAFETY DEVICE

It stops any inappropriate use when the carer is away.



## OTHER ACCESSORIES



### AC0838 IV POLE

The IV pole is provided with a supporting stand made of stainless steel which can be applied to any of the four corners of the upper structure of the shower trolley. This support can be adjusted as to its rotation. The IV pole can be adjusted as to its height. Pole dimensions: 110 x 23,5 cm. Support dimensions: 28 x 2 cm



### BUMPERS

All Chinesport shower trolleys are fitted with swivelling bumpers integrated into the deck structure.





Our Research & Development department has taken particular care in choosing high quality materials, components and special processes, ensuring the product is long lasting and rust-free. The video offers an overview on main features, details and possible adjustments.

# Shower trolleys for paediatric use

## Pos. 3-4 Intended use

- ☐ **3 1** Shower trolley 31  
for paediatric use

## Pos. 5 Frame height adjustment

- ☐ **C** Hydraulic, with pedal  
☐ **D** Electric, with hand control  
☐ **G** Electric, foot switch on its base

## Pos. 6 Mattress inside width

- ☐ **1** Inside width 60 cm

## Pos. 7 Mattress length extension

- ☐ **1** No length extension mattress  
☐ **2** 7 cm length extension mattress  
☐ **3** 14 cm length extension mattress

## Pos. 8 Backrest

- ☐ **1** Non-tilting  
☐ **2** Tilting

## Pos. 9 Wheels and braking system

- ☐ **1** 12.5 cm wheels diameter and individual brake on each wheel

L	D							
1	2	3	4	5	6	7	8	9





150 kg



Model shown - code LD31C1121

Dimensions are in cm

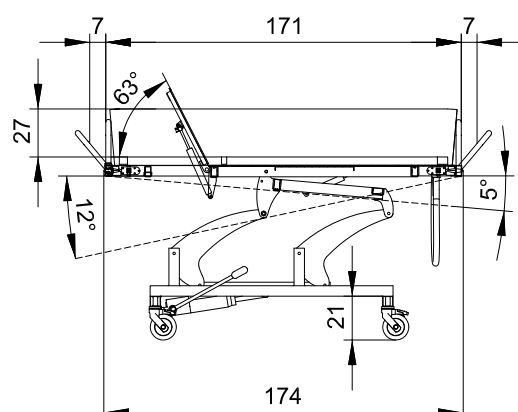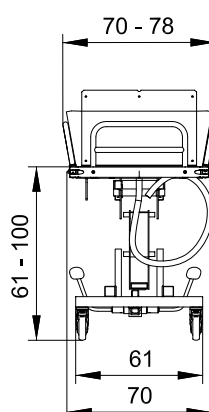

The backrest tilt can be adjusted up to a 63° angle via ergonomic control.

## TECHNICAL DATA

Height adjustment system	Hydraulic / Electrical with Battery
IP protection	X5
Supply	battery charger 110-220V 50/60 Hz 500mA
Electrical safety	Internal supply 24V
Safe working load	170 kg (body weight 150 kg)
Weight	80 kg
Wheel diameter	Ø 12.5 cm
Braking system	Individual on each wheel
Directional system	√
Side rails width reduction	√
Trendelenburg tilt	12°
Tilt for water discharge	5°
Discharge hose length	150 cm

## ACCESSORIES

### SHOWER MATTRESS ACCESSORIES

AC0889 SHOWER MATTRESS - S  
 AC1236 WEDGE PILLOW  
 AC1237 SEMI-CYLINDRICAL CUSHION  
 AC1126 SHOWER MATTRESS STRAPS

### ELECTRICAL ACCESSORIES

AC0071 SPARE BATTERY  
 AC0318 WALL BATTERY CHARGER  
 AC0578 SAFETY DEVICE

### OTHER ACCESSORIES

AC0838 IV POLE

## Pos. 3-4 Intended use

- ☐ **3 2** Shower trolley 32 for adults

## Pos. 5 Frame height adjustment

- ☐ **C** Hydraulic, with pedal  
☐ **D** Electric, with hand control  
☐ **G** Electric, foot switch on its base

## Pos. 6 Mattress inside width

- ☐ **2** Inside width 68 cm  
☐ **3** Inside width 80 cm

## Pos. 7 Mattress length extension

- ☐ **1** No length extension mattress  
☐ **2** 7 cm length extension mattress  
☐ **3** 14 cm length extension mattress

## Pos. 8 Backrest

- ☐ **1** Non-tilting  
☐ **2** Tilting

## Pos. 9 Wheels and braking system

- ☐ **1** 12.5 cm wheels diameter and individual brake on each wheel  
☐ **2** 15 cm wheels diameter and central brake



The rail width can be reduced to 78 cm to pass standard doors in a comfortable way.

L	D								
1	2	3	4	5	6	7	8	9	




150 kg

12



Model shown - code LD32C2122

Dimensions are in cm

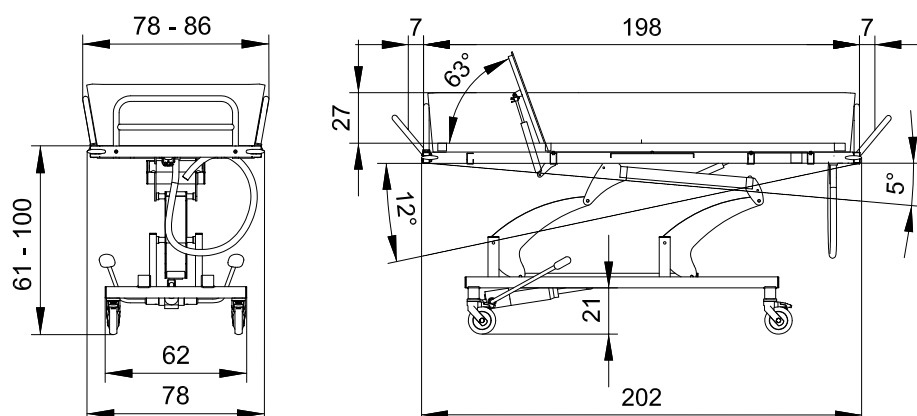

## TECHNICAL DATA

Height adjustment system	Hydraulic / Electrical with Battery
IP protection	X5
Supply	battery charger 110-220V 50/60 Hz 500mA
Electrical safety	Internal supply 24V
Safe working load	170 kg (body weight 150 kg)
Weight	80 kg
Wheel diameter	Ø 12.5 cm / Ø 15 cm
Braking system	Individual / centralized
Directional system	√
Side rails width reduction	√
Trendelenburg tilt	12°
Tilt for water discharge	5°
Discharge hose length	150 cm

## ACCESSORIES

### SHOWER MATTRESS ACCESSORIES

AC0890 SHOWER MATTRESS - M  
 AC0891 SHOWER MATTRESS - L  
 AC1236 WEDGE PILLOW  
 AC1237 SEMI-CYLINDRICAL CUSHION  
 AC1126 SHOWER MATTRESS STRAPS  
 AC0574 BACKREST EXTENSION

### ELECTRICAL ACCESSORIES

AC0071 SPARE BATTERY  
 AC0318 WALL BATTERY CHARGER  
 AC0578 SAFETY DEVICE

### OTHER ACCESSORIES

AC0838 IV POLE

# Shower trolleys *for bariatric use*

## Pos. 3-4 Intended use

- ☐ **3 3** Shower trolley 33  
for bariatric use

## Pos. 5 Frame height adjustment

- ☐ **D** Electric, with hand control

## Pos. 6 Mattress inside width

- ☐ **2** Inside width 68 cm  
☐ **3** Inside width 80 cm

## Pos. 7 Mattress length extension

- ☐ **1** No length extension mattress  
☐ **2** 7 cm length extension mattress  
☐ **3** 14 cm length extension mattress

## Pos. 8 Backrest

- ☐ **1** Non-tilting

## Pos. 9 Wheels and braking system

- ☐ **2** 15 cm wheels diameter and  
central brake

L	D								
1	2	3	4	5	6	7	8	9	





225 kg



Model shown - code LD33D2312

Dimensions are in cm

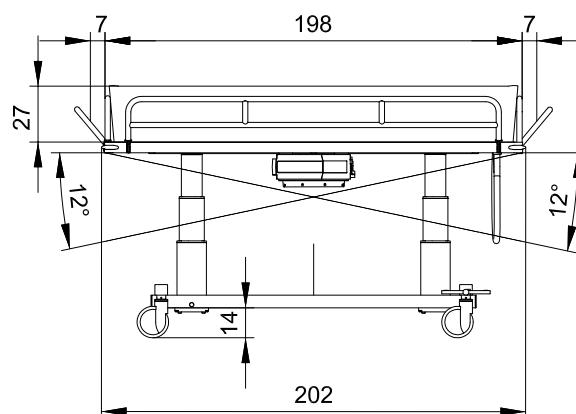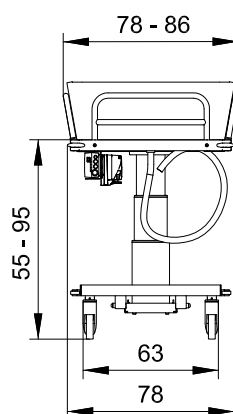

Rails can be dropped with one hand. Mattress height is 27 cm.

## TECHNICAL DATA

Height adjustment system	Electrical with Battery
IP protection	X5
Supply	battery charger 110-220V 50/60 Hz 500mA
Electrical safety	Internal supply 24V
Safe working load	245 kg (body weight 225 kg)
Weight	106 kg
Wheel diameter	Ø 15 cm
Braking system	Centralized
Directional system	√
Side rails width reduction	√
Trendelenburg tilt	12°
Tilt for water discharge	12°
Discharge hose length	150 cm

## ACCESSORIES

### SHOWER MATTRESS ACCESSORIES

AC0890 SHOWER MATTRESS - M

AC0891 SHOWER MATTRESS - L

AC1236 WEDGE PILLOW

AC1237 SEMI-CYLINDRICAL CUSHION

AC1126 SHOWER MATTRESS STRAPS

### ELECTRICAL ACCESSORIES

AC0071 SPARE BATTERY

AC0318 WALL BATTERY CHARGER

### OTHER ACCESSORIES

AC0838 IV POLE



## 01933 SHOWER MATTRESS

It allows bed-ridden, elderly or disabled patients to be washed and cared for in their personal hygiene with running water directly on their bed, without any need for transfer/lift. Inside the "Shower Mattress", the patient can be showered with water coming from the tap via an extension hose, or alternatively, using the special pressurized water tank provided with flexible hose and shower head. Once finished, the water is drained into a suitable container. Complete with plug and flexible hose for quick water drain.

Dimensions: 186 x 80 x 17h cm



12



1. Spread the shower mattress out on one side of the bed, next to the person being assisted.



2. Roll the person assisted on one side and slide the rigid side panel underneath him, using the suitable white handles if necessary.



3. Roll the person assisted onto their back and put the mattress in the centre of the bed, using the relevant handles and then hook up the corners.

## ACCESSORIES



## 01938 SHOWER EXTENSION HOSE

Flexible, 6 m hose complete with button-operated shower head and tap adapter.



## 01939 CARRY BAG

A container for holding and carrying the mattress.



## 01940 WATER TANK









Pressurized water tank complete with flexible hose and shower head operated by button.

## OTHER SPARE PARTS AND ACCESSORIES












NR	PICTURE	CODE	DESCRIPTION	LD31	LD32	LD33
1		<b>SP0019</b> WHEEL 125	Caster for base models code LD31 and LD32 in case of configuration with individual brakes. 125 mm diameter. (single unit)	•	•	
2		<b>SP0020</b> STEERING WHEEL 125	Caster with direction block for base models code LD31 and LD32 with individual brakes. 125 mm diameter. (single unit)	•	•	
3		<b>SP0017</b> WHEEL 150	Caster for base models code LD31 and LD32 with centralized brakes configuration. 150 mm diameter. (single unit)		•	•
4		<b>SP0018</b> STEERING WHEEL 150	Caster for base models code LD32 and LD33 in case of configuration with centralized brakes. 150 mm diameter. (single unit)		•	•
5		<b>SP0015</b> LEVER	Right and left pedal for models with centralized brake system. (a couple)		•	•
6		<b>SP0001</b> CONTROL UNIT	Control unit for shower trolley for bariatric use with two motors base code base LD33			•
7		<b>SP0030</b> CONTROL POWER LD33	Power unit for shower trolley base code LD33, with height and tilt adjustment functions to obtain the Trendelenburg position or for water drainage.			•
8		<b>SP0012</b> CONTROL POWER LD31-LD32	Power unit for shower trolley with one motor base code LD31 and LD32, with height adjustment function of the frame.	•	•	
9		<b>SP0014</b> MOTORIZED COLUMN	Motorized column for base model code LD33. Eventual substitution must be carried out only by qualified personnel.			•
10		<b>SP0013</b> ELECTRIC ACTUATOR	Electrical actuator for models with base code LD31 and LD32, for models with paediatric and adult use.	•	•	
11		<b>SP0016</b> HAND CONTROL 2	Hand control for models with one motor for adjusting the frame height and getting the Trendelenburg position, can be used also as a second control and therefore it includes a double cable system.			•

## OTHER SPARE PARTS AND ACCESSORIES

12

NR	PICTURE	CODE	DESCRIPTION	LD31	LD32	LD33
12		<b>AC0039</b> HAND CONTROL 1	Hand control for models with one motor for adjusting the frame height, can be used also as a second control and therefore it includes a double cable system.	•	•	
13		<b>AC0071</b> SPARE BATTERY	It is recommended in order to avoid interruptions during adjusting height whenever the main battery is being recharged.	•	•	•
14		<b>AC0083</b> FOOT SWITCH	Pedal control for models with one motor for frame height adjustment; it can also be used as a second control, it includes a cable.	•	•	
15		<b>AC0575</b> FOOT SWITCH HOLDER	Support designed for the use of the pedal control AC0083. It is useful in case of transferring the shower trolley and when transferring a patient.	•	•	
16		<b>SP0009</b> HYDRAULIC PISTON	For hydraulic adjustments of the frame height. Such version is considered for LD31 and LD32 models, for models with paediatric and adult use. (single unit)	•	•	
17		<b>SP0010</b> DOUBLE PEDAL	Right and left to be used with hydraulic piston SP0009. To be inserted in the right and left pedals of the hydraulic piston to ease the manual operation of pressure and release. (a couple)	•	•	
18		<b>SP0005</b> RUBBER COVERS	To be inserted in the right and left pedals of the hydraulic piston to ease the manual operation of pressure and release. (a couple)	•	•	
19		<b>SP0002</b> HANDLE	It is used to adjust the tilt of the shower trolley both for draining and for obtaining the Trendelenburg position. (single unit)	•	•	
20		<b>SP0011</b> GAS SPRING	The adjustments for draining and Trendelenburg are both activated by means of gas spring connected with a dedicated cable to the SP0002 handle. (single unit)	•	•	
21		<b>SP0026</b> CABLE	Connecting cable to activate the dedicated gas spring, which is used for adjusting tilt of the shower trolley in both directions. (single unit)	•	•	
22		<b>AC0892</b> DRAINING KIT	It includes a draining hose of 150 cm long and a drain to be applied to the mattress; (if needed a different length of the draining hose can be requested)	•	•	•


## OTHER SPARE PARTS AND ACCESSORIES

NR	PICTURE	CODE	DESCRIPTION	LD31	LD32	LD33
23		<b>SP0041</b> MATTRESS SMALL	Mattress to be used for the shower trolley paediatric model base code LD 31, 160 x 60 cm size. It does not have a draining hose or plug which have to be requested separately.	•		
24		<b>SP0035</b> MATTRESS MEDIUM	Mattress to be used for the shower trolley for adult and bariatric models base code LD32 and LD33, 188 x 68 cm size. It does not have a draining hose or plug which have to be requested separately.		•	•
25		<b>SP0042</b> MATTRESS LARGE	Mattress to be used for the shower trolley for adult and bariatric models base code LD32 and LD33, 188 x 80 cm size. It does not have a draining hose or plug which have to be requested separately.		•	•
26		<b>SP0021</b> PLASTIC COVER CAP 50X50	Set of caps 50 x 50 mm used for models with individual brake system. To be applied onto the base frame. It is for use in one trolley only. (8 units per set)	•	•	
27		<b>SP0022</b> PLASTIC COVER CAP 60X50	Set of caps 60 x 50 mm used for models with individual brake system. To be applied onto the base frame. It is for use in one trolley only. (4 units per set)		•	•
28		<b>SP0023</b> PLASTIC COVER CAP 60X50 WITH HOLE	Set of caps 60 x 50 mm with a 17mm Ø hole used for models with centralized brake system. To be applied onto the base frame. It is for use in one trolley only. (4 units per set)		•	•
29		<b>SP0024</b> PLASTIC COVER CAP 50X30	Set of caps 50 x 30 mm. To be applied onto the central part of the frame. It is for use in one trolley only. (8 units per set)	•	•	
30		<b>SP0025</b> PLASTIC COVER CAP 80X30	Set of caps 80 x 30 mm. To be applied onto the central part of the frame, half-moon levers. It is for use in one trolley only. (2 units per set)	•	•	
31		<b>SP0027</b> HEXAGON NUT COVER M12	Set of nut covers M12 used in the central part of the frame and the lever components. It is for use in one trolley only. (20 units per set)	•	•	
32		<b>SP0028</b> HEXAGON NUT COVER M8	Set of nut covers M8 used in the upper part of the frame. It is for use in two trolleys. (20 units per set)	•	•	
33		<b>SP0029</b> HEXAGON NUT COVER M10	Set of nut covers M10 used in the hydraulic piston components. It is for use in one trolley only. (4 units per set)	•	•	



## LD31 SHOWER TROLLEY 31 for paediatric use


12

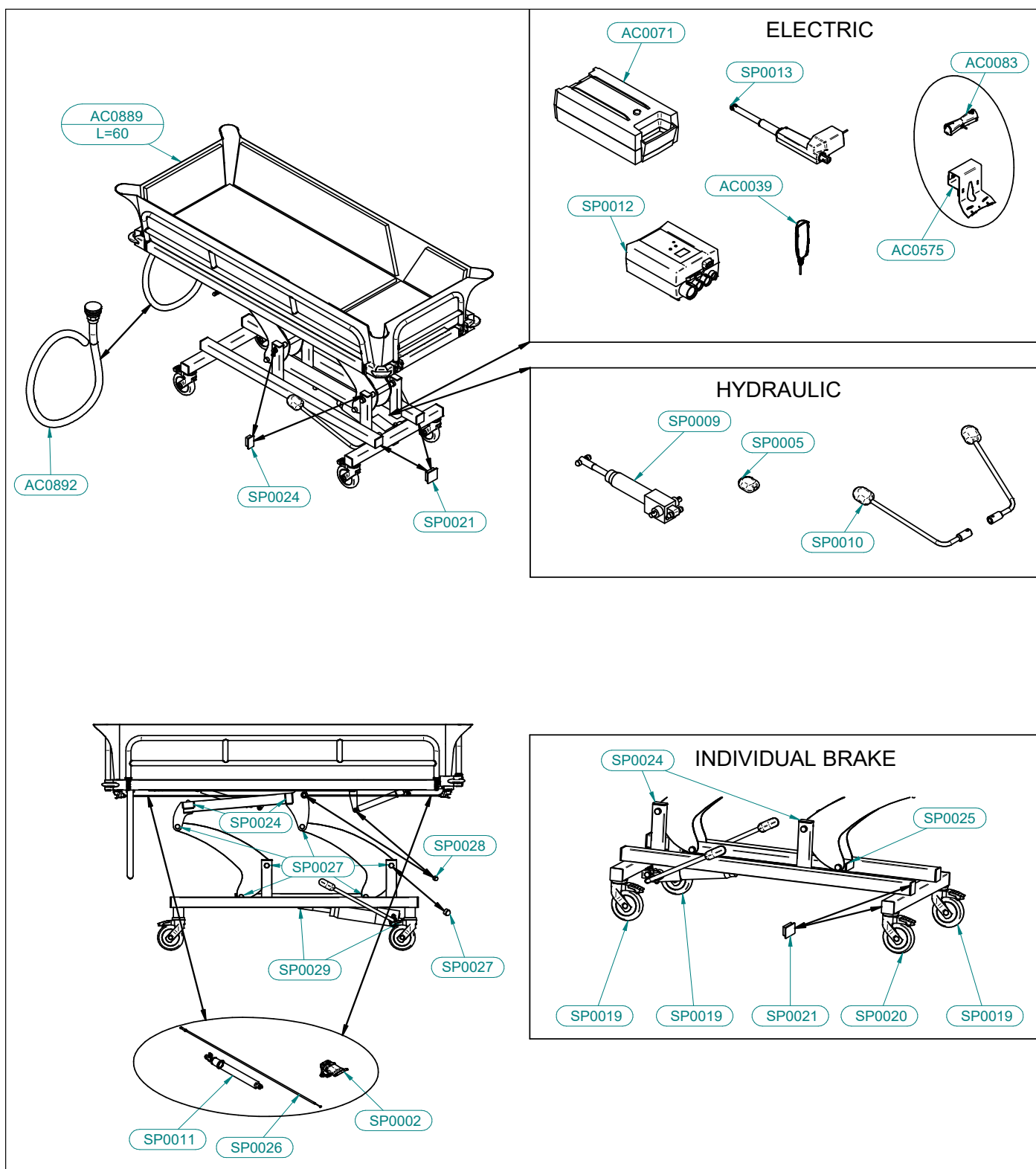

# Shower Trolleys *Exploded views*


## LD32 SHOWER TROLLEY 32

for adults

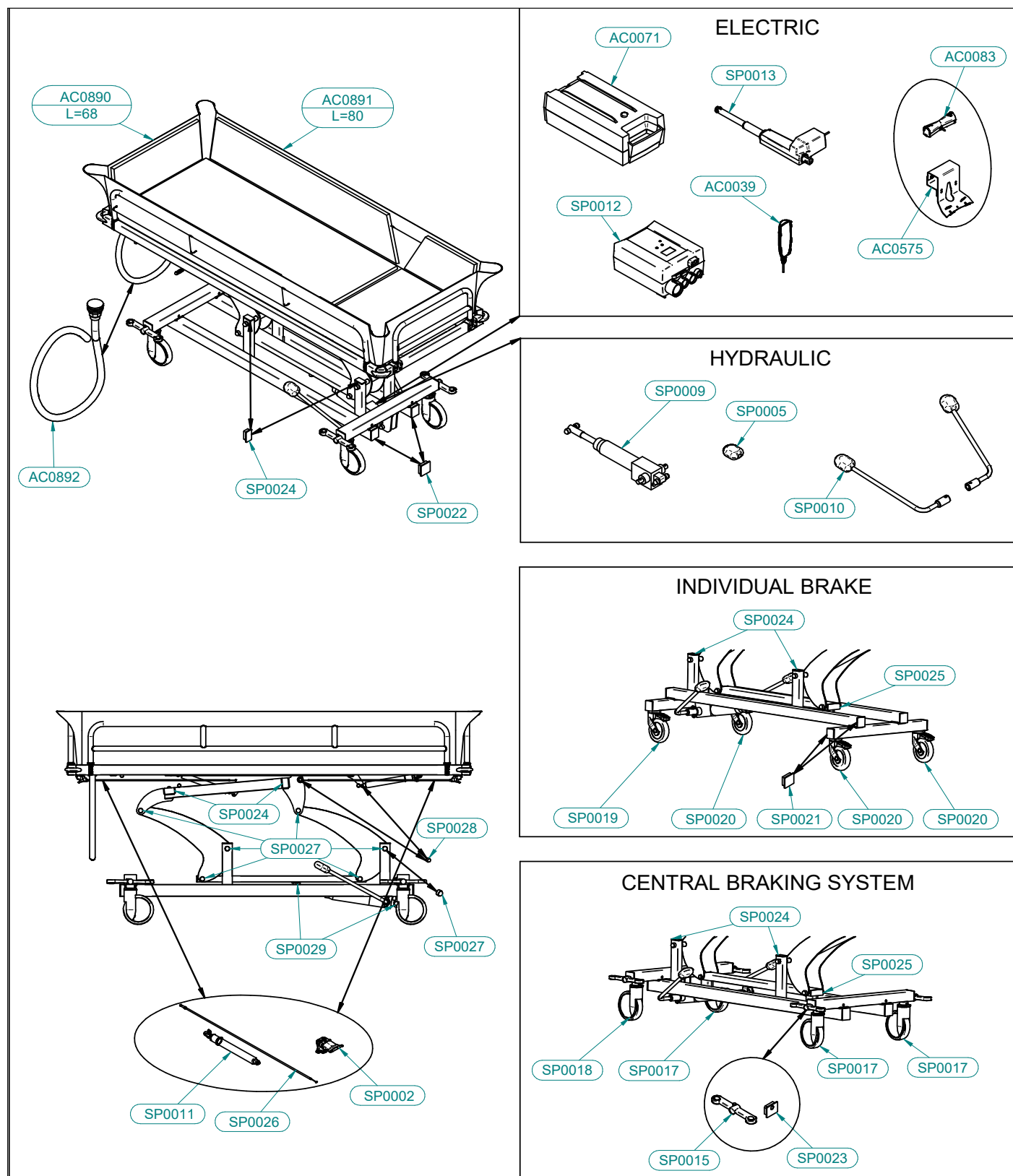

## LD33 SHOWER TROLLEY 33 for bariatric use

12

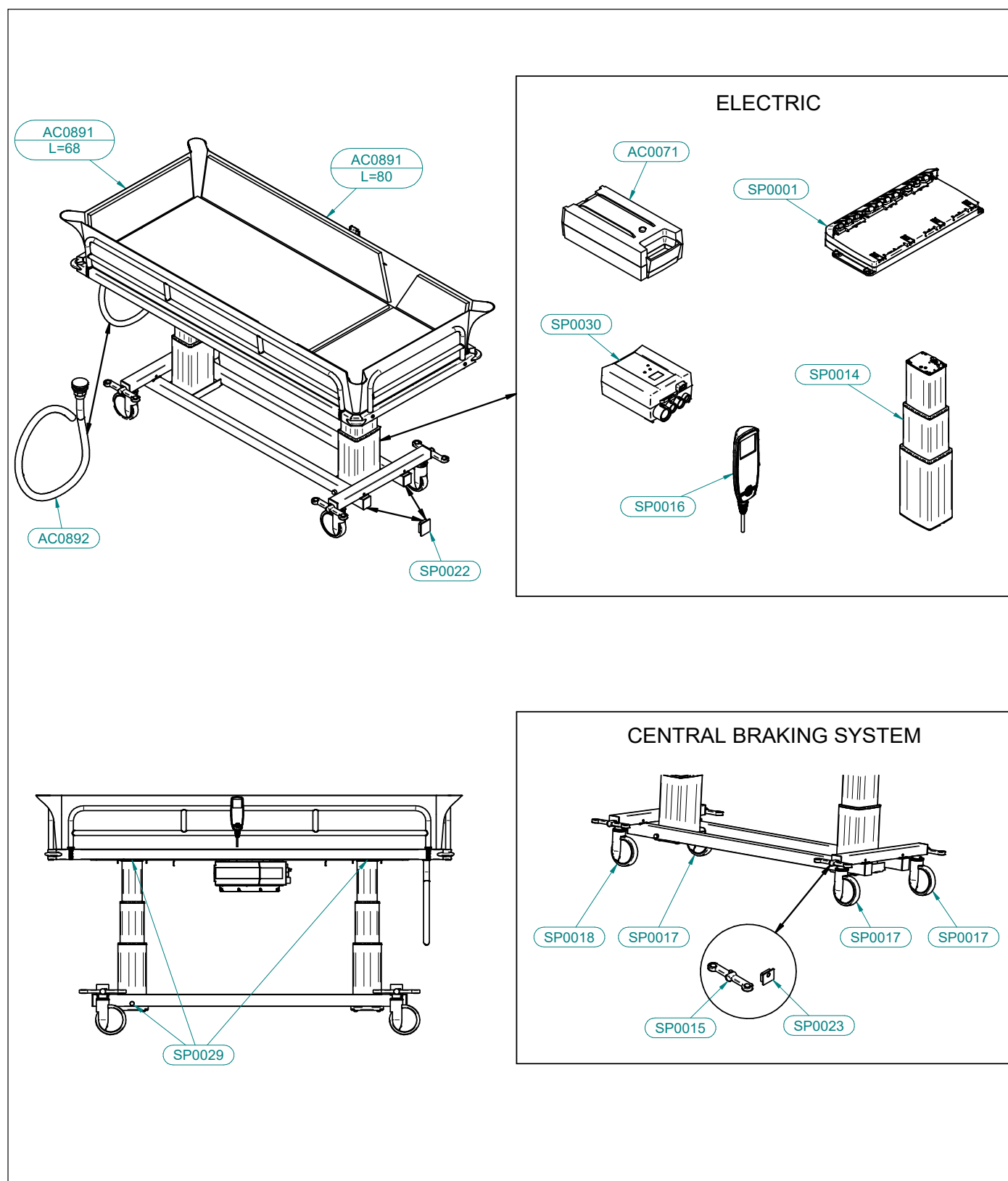



Chinesport is based in Udine, Italy, between the Alps and Venice. For over 40 years we have been dedicated to healthy posture for healthy movement. The root of our company name refers to the Italian word “chinesiterapia”, or movement therapy. We strongly believe and adhere to “movement culture” as a way to prevent and cure injury and disease.

Today we are a global leader in developing and manufacturing rehabilitation equipment and assistive devices. We have excellent and long-standing business relationships in almost 80 countries worldwide. The Chinesport general product catalogue contains over 1.000 innovative, high-quality products. New catalogue editions that include the latest product innovations and trends are regularly published. Our own medical-scientific training and educational program is continuously expanding and caters for all specialised rehabilitation fields. As an organisation, we have been working with a certified quality management system and in compliance with international ISO 9001 and ISO 13485 standards since 1998.



CHINESPORT spa - Via Croazia, 2 - 33100 Udine - Italy  
Phone 0432 621 621 - Fax 0432 621 620 - [export@chinesport.it](mailto:export@chinesport.it)

[www.chinesport.com](http://www.chinesport.com)